



# HL-SK-300

## Sliding Stretcher Cart

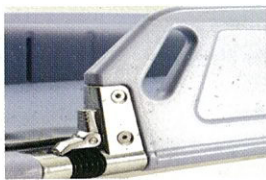

### Side Rail Holder

Safety on patients get higher as the double locking system has been applied

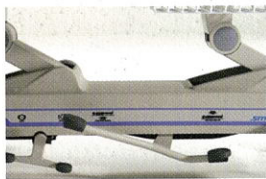

### Pedal

Center-concentrated brake pedal and height-control pedal



### O2 Tank hanger

(Option)



### Sliding Mattress

Pull down the side rail after equalizing the height of the car with that of the bed. Then, patients can be easily transferred as the mattress moves towards the bed. Plus, a clutch is installed inside the steering wheels, so that steering wheels idle when the mattress finishes its movement. Like this, it can be prevented from being broken down.

### Sliding Stretcher Cart

Model	Size (LxDxH / mm)	Remark
HL-SK-300	2020x660x530~830	

### SPEC.

- Frame : Steel pipe
- Body : Aluminum
- Pump : Oil hydraulic pump
- Finish : Electrostatic powder coating(Use of the epoxy-polyester)
- Casters : 6inch swivel caster(Central locking system)
- Top : P.P
- Base : Steel pipe
- Mattress : Sliding mattress

### Std. Acc.

- I.V. pole 1EA

