








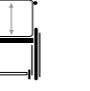


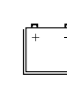






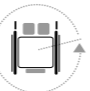



Technical information Navix

|  |  |  |  |  |  |  |  |  |  |  |
|--|---|---|---|---|---|---|--|---|---|---|
| total width | seat width | seat inclination | back inclination | height armrests | seat height | total height | total length | seat depth | back height | max. users weight |
| 570 | 400 | 5°/7,5°/10° | 10°/20°/30°/40° | 260-350 | 400/450/500 | 995 | 945 | 380-500 | 580 | 130 |
| 570 | 450 | 5°/7,5°/10° | 10°/20°/30°/40° | 260-350 | 400/450/500 | 995 | 945 | 380-500 | 580 | 130 |
| 620 | 500 | 5°/7,5°/10° | 10°/20°/30°/40° | 260-350 | 400/450/500 | 995 | 945 | 380-500 | 580 | 130 |

Complies to crash test ISO 7176-19

Deviations in sizes and weight (1,5 mm and 1,5 kg) are possible

|  |  |  |  |  |  |  |  |  |  |
|--|---|---|---|---|---|---|---|---|---|
| weight | batteries | power | speed | driving range | max. slope upwards | max. slope downwards | ground clearance | turning radius | wheels |
| 79 | 2x12V 38Ah | 2 x 220 W | 6 km/h | 20 km | 9° | 9° | 60 | 1250 | front: air 317,5x57,15 rear: PU 200x50 |

Class A wheelchair for indoor use. - Not intended as passenger seat during transport.

Options

- Plexi Table
- LED Lights
- Anti-tipping
- Inclinable footplate
- Sideward controller
- Different seat cushions
- Different back cushions
- Side support
- Headrest
- Electrical legrests
- Controller adjustable in height

Extra option: LIFT

300 mm lift for a better SOCIAL INTEGRATION

The lift allows me to interact with people around me, enhancing my mobility and social integration.

300 mm lift for a better ACCESSIBILITY

- Kitchen cupboards
- Cash dispensers
- Bar
- etc...



Colors



* Additional colors on demand

N.V. Vermeiren N.V.

Vermeirenplein 1-15
B-2920 Kalmthout
Belgium

Tel.: +32 (0) 3 620 20 20
Fax: +32 (0) 3 666 48 94

www.vermeiren.com - info@vermeiren.be



VERMEIREN
Opening up new horizons in mobility

For more information, contact your local dealer



Subject to technical alterations



Better maneuvering with a compact front-wheel drive



- Better overcoming obstacles
- Better taking sharp turns
- Easy to transfer
- Compact - fits in a small car
- Adjustable to your needs



model shown is for illustration only



Opening up new horizons in mobility

Information for the visually impaired is available upon request. - R.E.: N.V. Vermeiren N.V., Vermeirenplein 1/15 - 2920 Kalmthout - Belgium - 12/2011



Better maneuvering with a compact front-wheel drive

The Navix allows for sharper turning indoor and for better overcoming obstacles like tresholds when entering the local grocery store.

The user friendly Navix is adjustable to individual needs and allows for easy transfers to bed, bath, toilet, car etc...

The compact Navix fits in a small car for even more mobility.



01

Better maneuvering with front-wheel drive

- sharper turning: front-wheel drive vs rear-wheel drive
- better overcoming obstacles

“Maneuvering inside my house was never so easy with my Navix.”

Compact

- easy to transport in small car
- easy to reassemble
- total width for Navix limited to 570 mm.

02



03

User friendly

- reclinable backrest and seat
- adjustable seat depth
- total width indication
- easy to transfer (armrest)
- adjustable seat height/angle
- adjustable armrests in height and width
- adjustable depth controller (optional in height and fold sideways)



04

Outdoor usage

- ground clearance of 60 mm
- better grips on outdoor surfaces with profiled large front wheels

