



CGM Hospital Furniture

3. Bedside Cabinet / Mobile IV Stand / Overbed Table

1. Bedside Cabinet


- Clean lines in environmental friendly material.
- Firm, anti-rust, damp-proof, heat-resistant, scratch-resisting and anti-acidity.

- Smooth table with fence around to avoid things falling off and easy to clean.
- Available colors of Maple and Cherry in wooden grain for choices.
- DIY assembling or complete-assembled available for requests.

Maple Color	Specification Cherry Color	Description
		<ul style="list-style-type: none"> ● Dimension : L55 x W48x H77 cm ● Built-in both side towel hangers ● With shoes rack (Reticular design for good ventilation) ● Drawer, closet and sliding board available ● Caster : 4 x 2" caster (2 with brakes)


2. Mobile IV Stand

- Outstanding features in design, appearance and material to upgrade patients' security and comfort.
- Stand base in triangular shape enables patients walk comfortably and safely.
- PU foam base avoids unexpected collision for patients when walking.
- Stand base set in weight is capable of descending the accident of turnover.

Specification	Description
	<ul style="list-style-type: none"> ● Flexible Height adjustable from 160 to 200cm ● Material : <ul style="list-style-type: none"> - Telescoping tube, IV hookers and pushing handle in stainless (SUS 304) - Stable stand base and pulling handle in PU foam ● Screw knob for height adjustment, slow descending via inside buffer mechanism ● IV hanging hooks (2 or 4) available for choice ● Caster : 1 x 4" steering and 2 x 4" swivel

3. Overbed Table

- Distinguished design in the aim of stable, durable, reliable and be a functional partner for patients.
- The widening lifting column with gas-spring assisted perfectly presents its easy-operation and stable functions.
- Table in the material of anti-rust, damp-proof, heat-resistant, scratch-resistant and even anti-acidity.

Specification	Description
	<ul style="list-style-type: none"> ● Dimension : <ul style="list-style-type: none"> - Table : L90 x W 42cm - Adjustable Height : 76 ~ 105cm ● Safety Weight Capacity : 30Kg ● Material : <ul style="list-style-type: none"> - Table : Chip Board with laminate in wooden grain around with PU rim. - Lifting Column : Aluminum, Steel - Base : Aluminum with Powder Coating ● Caster : 2 x 1" and 2 x 2"

CGM reserves the right to make changes without notice in design and specification.

CGM Chang Gung Medical Technology Co., Ltd.

Headquarter: 11F., No.201-4, Dun-Hwa North Rd., Taipei 105, Taiwan, R.O.C.
Tel: 886-2-87122948 Fax: 886-2-25140620

E-mail: cgmcc@gmcc.com.tw http://www.cgmcc.com.tw













CGECF1012-1



Chang Gung Medical Technology Co., Ltd.

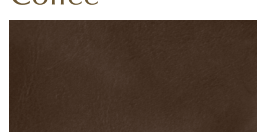

1. Care Recliner

Specification		
Safety Knob  The knobs situated at both sides of the recliner. Slightly pulling either one of the knobs can transform the recliner to a bed or a chair easily.	Gas Spring  There is a Gas Spring set at the back of the recliner. This supporting device can effectively assist users to take less energy when transforming the recliner.	Total Smooth Bed Surface  Ergonomic sleeping surface with rotating side arms provides comfortable and bigger space to rest.
Operating Handle  Situated at the both sides of the recliner. Easy to move or convert the recliner.	Braking Caster  80mm PU type braking caster. Step on the brake switch at the top of the caster can brake/un-brake the caster.	Side Arms  Foam side arms provide a comfortable touch. There are 3 sizes available for different requirements of wards.
Overall Size		
Large Chair: 90L x 74W x 107H cm Bed: 217L x 74W x 45H cm Sleeping surface: 202L x 60W cm	Medium Chair: 87L x 69W x 100H cm Bed: 200L x 69W x 45H cm Sleeping surface: 190L x 55W cm	Small Chair: 86L x 60W x 98H cm Bed: 196L x 60W x 45H cm Sleeping surface: 185L x 52W cm
Color		
Coffee 	Beige 	Gray 
Pearl Red 		
* Please note that color choices only base in a considerable quantity.		

2. Medical Air Unit / Bed Locator

1. Medical Air Unit

- Each medical unit is equipped with three sets of Oxygen, Medical Air and Vacuum outlets services and sufficient rooms for chest bottle, Ambu, suction tubes, swabs and other medication. The medical air unit helps to provide patients with better care. The sliding door is easy for manual use and can effectively keep goods clean when sliding door is close.
- The elegant exterior design can be applied when designing a pleasant V.I.P. patient room.




White Color	Specification Wooden Color	Description
		<ul style="list-style-type: none"> ● Material: Special ABS container and aluminum outlets retainer ● Dimension: W528 x H657 x T182 mm ● Optional Accessories: Ambu Holder, Chest Bottle Holder * The unit does not include Oxygen, Medical Air, and Vacuum outlet.

2. Bed Locator

- The bed locator system not only facilitates positioning but also protects walls and equipments from unexpected collision during transportation.
- The Locator is capable of setting up with different devices, such as switches, plugs, electrical communication, night light and nurse call on both sides.

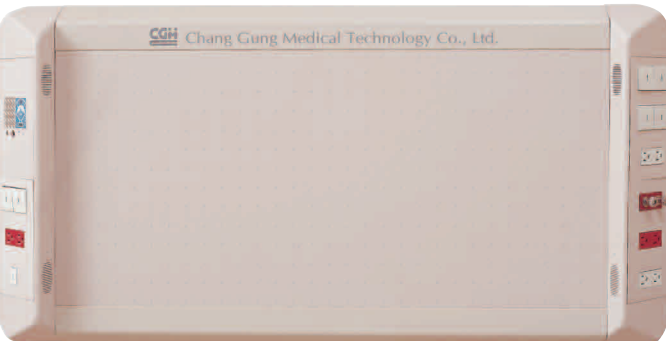

Specification
 <ul style="list-style-type: none"> ● Material: ABS frame and metal interior structure ● Dimension: W1208 x H612 x T65 mm * The system does not include switch, receptacles, nurse call, and communications connectors.