



#### 4. Various Optional Accessories

Small things bring more convenient caring.

Chang Gung Clinical Cart is outstanding not only in various ABS modulation assembling possibility, but also always keeps up with world-trend technology to meet users' requirements worldwide. Some optional accessories we provide as below:

Working Tray in Different Design			
Front Working Tray (Pulling Type)	Side Working Tray (Pulling Type)	Side Working Tray (Lifting Type)	
Electronic Lock & Key Lock			
Electronic Lock (Bracelet)	Electronic Lock (Card)	Key Lock	
ABS Drawers with ABS & Aluminum Tracks			
DIY Dividable Drawer	ABS & Aluminum Track		
More Accessories			
Cardiac CPR Massage Board	Suction Tube Storage	Foot-Operated Waste Bin	Needle Container
Four Medicine Bottle Dispensers			
A Type	B Type	C Type	D Type
Monodose Cassette Sizes			
Small	Medium	Large	
<ul style="list-style-type: none"> <li>Total 10 cassettes for a layer</li> <li>Size : 18L x 10.5W x 7H cm</li> </ul>	<ul style="list-style-type: none"> <li>Total 8 cassettes for a layer</li> <li>Size : 22L x 12W x 7H cm</li> </ul>	<ul style="list-style-type: none"> <li>Total 6 cassettes for a layer</li> <li>Size : 22L x 16W x 7H cm</li> </ul>	

Specifications	
	<b>1 Hanging Basin</b> <ul style="list-style-type: none"> <li>Size : 45L x 14.7W x 12H cm</li> <li>Available to place cans and bottles.</li> <li>Adjustable dividers available.</li> </ul>
	<b>2 ABS Drawer (Small)</b> Size : 52L x 33W x 8H cm Note : The install ABS tray can be taken out for cleaning and sterilizing.
	<b>3 ABS Drawer (Medium)</b> Size : 52L x 33W x 16H cm
	<b>4 ABS Drawer (Large)</b> Size : 52L x 33W x 24H cm
	<b>5 Hanging Wirebasket</b> Size : 38.7L x 14W x 16H cm Note : Available to place cotton tissue and operating gloves.
	<b>6 Bottom Basket</b>
	<b>7 Charger and Battery Set</b>
	<b>8 Foot-Operated Waste Bin</b> Size : 40L x 17W x 47.8H cm Note : Enable to place infected and non-infected waste.
	<b>9 Needle Container</b> Size : 18L x 14.5W x 18H cm(4.5liter) 18L x 7.3W x 18H cm(2.2liter)
	<b>10 Medicine Bottle Dispenser</b>
	<b>11 Cardiac CPR Massage Board</b> Size : 57L x 45W x 0.8T cm
	<b>12 PC Swing Arm</b> The arm is workable to different kinds of computers, e.g Panel PC or Notebook.
<b>Dimension</b> <ul style="list-style-type: none"> <li>Body : 75L x 60W x 100H cm (standard)</li> <li>Top Surface : 60L x 45W cm</li> <li>Caster : 5" x 4PC (2 with brakes)</li> </ul>	
<b>Exploded View</b> <p>A flexible assembling design to meet different requests.</p>	

CGM reserves the right to make changes without notice in design and specification.

**CGM Chang Gung Medical Technology Co., Ltd.**

Headquarter: 11F., No.201-4, Dun-Hwa North Rd., Taipei 105, Taiwan, R.O.C.  
 Tel: 886-2-87122948 Fax: 886-2-25140620

E-mail: cgmc@cgmc.com.tw http://www.cgmc.com.tw



## CGM Clinical Cart



Chang Gung Medical Technology Co., Ltd.



### 1. Configurable Design for Different Caring Need

Features about Clinical Cart...


- ABS Plastic + Aluminum = Rust-Proof, Quieter and Easy to Clean.**
  - ABS plastic and aluminum columns are rust-proof, easy-cleaning and maintenance.
  - Providing a quieter and peaceful environment.
- Multi-Functional Modulation = Maximise the Assembling Possibility.**
  - Aim to meet perfectly different caregiver's requirement by unlimited assembling methods with "Main Body" and "Optional Accessories"
  - "Main Body", all made by ABS, such as Drawers/Dose drawers in different sizes, Cardiac CPR massage board, Working Tray ;
  - "Optional Accessories", such as Hanging Basin, Hanging Wirebasket, Needle Container, IV Pole Holder, Suction Tube Storage and etc..
  - The height of the cart can be modified according to the medication needs.
- Self-Dividable ABS Drawer = Outstanding & Flexible Content Management.**
  - The inside of the drawer could be flexibly self-divided based on different nursing requirements.
  - The track of each drawer is all made by ABS material and aluminum. Durable use, rust-proof and avoid shape-changed.



### 2. Good Solution for Non-paper, Non-film & Electronic Medication Environment

Features about e-Clinical Cart...

- The Functional PC Swing Arm**
  - The PC Swing Arm could be set with different types of PCs, either with Panel PC or Notebook based on nursing requirement.
  - The main specialties of the PC Swing Arm are easy to adjust any required angle, position and height of the PC for different height and needs of caregivers.
- Long Usage, Stable & Uninterrupted Power Supply**
  - The Power Supply system in the e-Clinical Cart features in "Long Usage", "Stable Power" and "Uninterrupted System".
  - Aim to attend safe use, reliable, convenient and efficient electronic medication environment.
- Distinguished Optional Accessories**
  - Optional socket, power charge and power storage system complete the e-Clinical Cart function.
  - All of optional accessories are modulation-oriented to go with the cart perfectly.

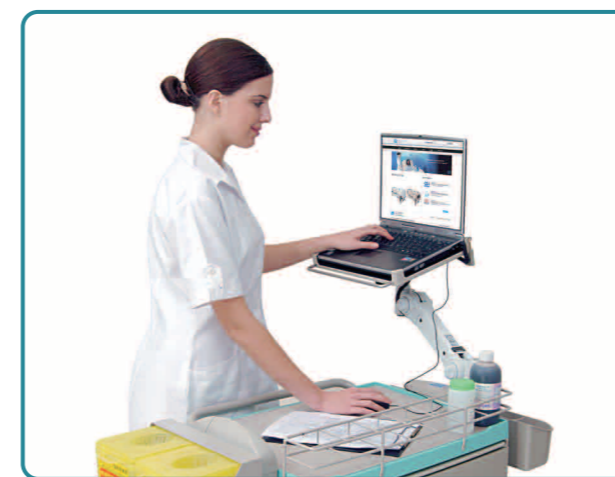

### 3. Multi-Functional & Assembling Cart

CGM carts aim to assemble as many configurations as you may require, such as e-Clinical Cart, Clinical Cart, UD Cart, Emergency Cart, Medication Cart and Layer Cart, etc. Some models showed here... (Please note that all our carts are able to be set up with e-Clinical Cart components.)

Model		
e-Clinical Cart Model A	e-Clinical Cart Model B	e-Clinical Cart Model C
Clinical Cart Model A	Clinical Cart Model B	2-Layer Cart Model A
UD Cart Model A	UD Cart Model B	UD Cart Model C
Emergency Cart Model A	Emergency Cart Model B	Emergency Cart Model C
Medication Cart Model A	Medication Cart Model B	Medication Cart Model C

Main Optional Accessories for e-Clinical Cart

**Power Display (Optional)**  
 Power Display is set on the top of Aluminum Column and shows three kinds of indicators: Power Condition, Low Power, and Recharging. To enhance better understanding on the power supply condition.

**Power Socket and Wire Collector (Optional)**  
 To provide two sets of triple hole socket, and Wire Collector with 2.5 meters long for power supply.

**Charger (Optional)**  
 Battery charger with "Red" means "recharging"; "Green" means "full".

**Battery Set (Optional)**  
 To provide different voltage of DC 19V / DC 12V for the choices.