







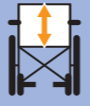



## Technical information Rapido

 total width	 folded width	 seat width	 height of armrests	 seat height	 total height	 total length	 seat depth	 back height	 weight (+batteries)
59	36	41	Min 20	53	98	105	44	51	38(+32)
64	36	46	ε	53	98	105	44	51	38(+32)
70	36	52	Max 25	53	98	105	44	51	38(+32)

Deviations in sizes and weight (± 1,5 cm and ± 1,5 kg) are possible

## Your guarantee of quality

Regarding quality, the Rapido complies with the **highest norms laid down by certified test institutes** in most European countries.

## The wheelchair must be adapted to you, not vice versa

The wheelchair has been developed as a **modern aid which can be individually adapted**. Aesthetic appearance and practical use go hand in hand. The wheelchair manufacturer Vermeiren has played an important role in this development.

Our basic principle is: the wheelchair must be adapted to you, not vice versa. This is why our wheelchairs have a **modular design**. All types of wheelchairs can be adapted, have numerous accessories and all options are **mutually interchangeable**. This is the result of our continuous research and more than 40 years of experience.

For more information, contact your local dealer

**N.V. Vermeiren N.V.**  
Vermeirenplein 1-15  
B-2920 Kalmthout  
België

Tel.: 03-620 20 20 • Fax: 03-666 48 94  
Also visit our web-site: [www.vermeiren.com](http://www.vermeiren.com)



# VERMEIREN

Opening up new horizons in mobility

Subject to technical alterations

# Rapido

Fast, Comfortable and  
fully adjustable



The Rapido combines **comfort and technology**. The adjustable back and seat will make the Rapido a part of a very personal experience. Up to date and easy to carry along due to a simple and fast foldable system.

R.E.: N.V. Vermeiren N.V., Vermeirenplein 1/15 - 2920 Kalmthout - Belgium - 10/2002



# Rapido



## Rapido briefly sketched

The Rapido can be best presented as an **all-in-one**. Comfort, technology and appearance go hand in hand. It is proven here once again that a marvellous piece of technology does not have to be a mastodon. Its 38kg makes the Rapido **highly manageable**.

The **anatomical seat and back** provide optimum support for the body and guarantee very high seating comfort for you. In addition, the armrests can be adjusted in both height and width or completely removed.

The **anti-tipping system** prevents the wheelchair from tipping over backwards.



## Comfort at your fingertips

The **control unit is collapsable** and allows you to move closer to the table. **Adjustable speeds** and a short braking distance will allow full control.

The **inclination of the back** (30°) and the seat (10°) provides optimal comfort in every environment. The footrests turn inward and outward or can be removed completely. Even the footplate can be adjusted in depth.



- 2x12V/38ah
- Average driving distance 30km
- Topspeed: max. 6km/u



## Flexibility with hidden talents

The Rapido engines are as **powerful** as its predecessors. Rechargeable, **maintenance-free** and a quick-action coupling ensuring that both batteries can be removed easily. Lifting a wheelchair in your car has never been easier.

A sensitive and easy manageable control unit ensures your **manoeuvre possibilities in even the most extreme circumstances**.



## Choice of 2 color combinations

Frame	Upholstery
Sparkling black C25	Black nylon
Metallic blue C27	Grey skai

## Options

- **Table**
  - Table control unit
  - Chin control
- **Lights and direction indicators**
  - Companion control
  - Anatomical seat and back
  - Adjustable legrests

