










Technical information Altitude

 total width	 seat width	 height armrests*	 seat height*	 total height	 total length	 seat depth	 back height*	 weight
67	38 40 42 44 48	23 25 27	45/47/49 51/52/53 54/55/56	119 121,5 124	112	46	75 77 79	25

Deviations in sizes and weight (1,5 cm and 1,5 kg) are possible

Construction

The Altitude is standard equipped with a 8 cm thick seating pillow (backheight : 77 cm). The height of the back and armrests are dependable on the thickness of the seating pillow (6-8-10 cm)(*).

The steel structure with a diameter of 32 mm and epoxy coating result in a very robust resistance against chemical and mechanical forces. The 5 wheels of 125 mm will give you a optimal steering comfort in almost any environment. The fifth steeringwheel will make it easy to steer while driving .

Seating pillow and back, made of foam rubber, will give additional comfort to the patient. The upholstery in Dartex is ideal in cases of extended sweating.

The wheelchair must be adapted to you, not vice versa

Our basic principle is: the wheelchair must be adapted to you, not vica versa. This is why our wheelchairs have a modular design.

All types of wheelchairs can be adapted, have numerous accessories and all options are mutually interchangeable. This is the result of our continuous research and more than 45 years of experience.

For more information, contact your local dealer

N.V. Vermeiren N.V.

Vermeirenplein 1-15
B-2920 Kalmthout
Belgium

Tel. : +32 (0)3 620 20 20
Fax: +32 (0)3 666 48 94

Visit our website: www.vermeiren.com



VERMEIREN

Opening up new horizons in mobility

Subject to technical alterations

Altitude

SAFETY & COMFORT



R.E.: N.V. Vermeiren N.V., Vermeirenplein 1/15 - 2920 Kalmthout - België - 09/2005

Options:




Footrest mounted with a spring which allow easy transfers in and out the wheelchair.

(Footrests are removable)



The pillow of the legsupport with side support and a separate spring, independent from the seat adjustment.


A simply removable table, adjustable in depth and mounted on a system which allows easy dismantling. The table will always be in a horizontal position no matter what the angle of the seat is.

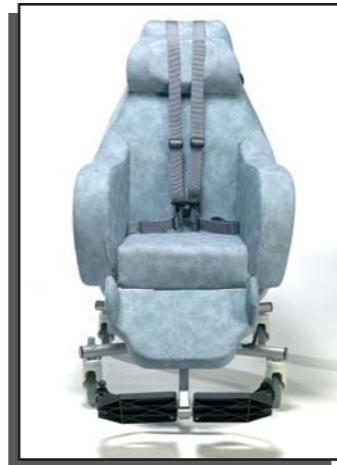

Safetybelt

Altitude

-4° to +36°

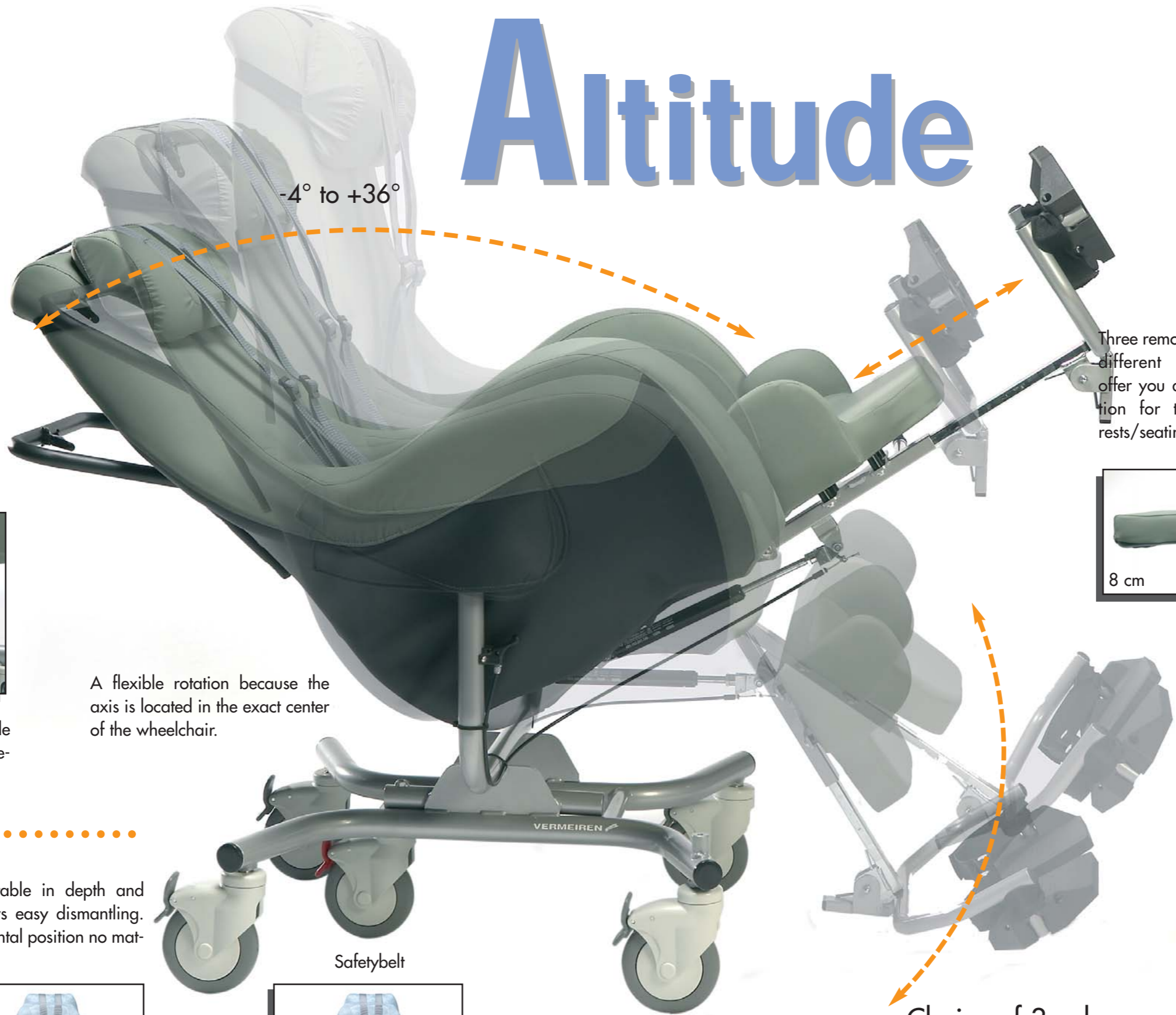

A flexible rotation because the axis is located in the exact center of the wheelchair.



Mounted with 4 wheels and one steering wheel.

Three removable pillows in different heights which offer you a height correction for the back/arm-rests/seating height.



6 cm



8 cm



10 cm



A spring allowing 3 positions creating a variable seating height from 45 to 56 cm.

Choice of 3 colours:



DARTEX Green



DARTEX Orange



SKAI Blue